



IHE Radiation Oncology

Test Tools Scenario Descriptions

Index


Index	2
Test Scenario 1 Archive actor	3
Test Scenario 1 Contourer actor	4
Test Scenario 2 Contourer actor	5
Test Scenario 3 Contourer actor	6
Equidistant Scenario 1 Contourer actor	7
Test Scenario 1 Dose Displayer	8
Test Scenario 2 Dose Displayer	9
Equidistant Test Scenario 1 Dose Displayer	10
Test Scenario 1 Dosimetric Planner	11
Test Scenario 2 Dosimetric Planner	12
Equidistant Test Scenario 1 Dosimetric Planner	13
Test Scenario 1 Geometric Planner	14
Test Scenario 2 Geometric Planner	15
Equidistant Test Scenario 1 Geometric Planner	16

Test Scenario 1 Archive actor

Description:

This test scenario provides the archive actor with all data objects and receives the data objects.

- Storage of the reference test data objects to the archive
- Retrieval of the reference test data object from the archive
- Comparison of the test data objects received, with the original reference data.
- Checks ensuring level 2 SCP conformance.

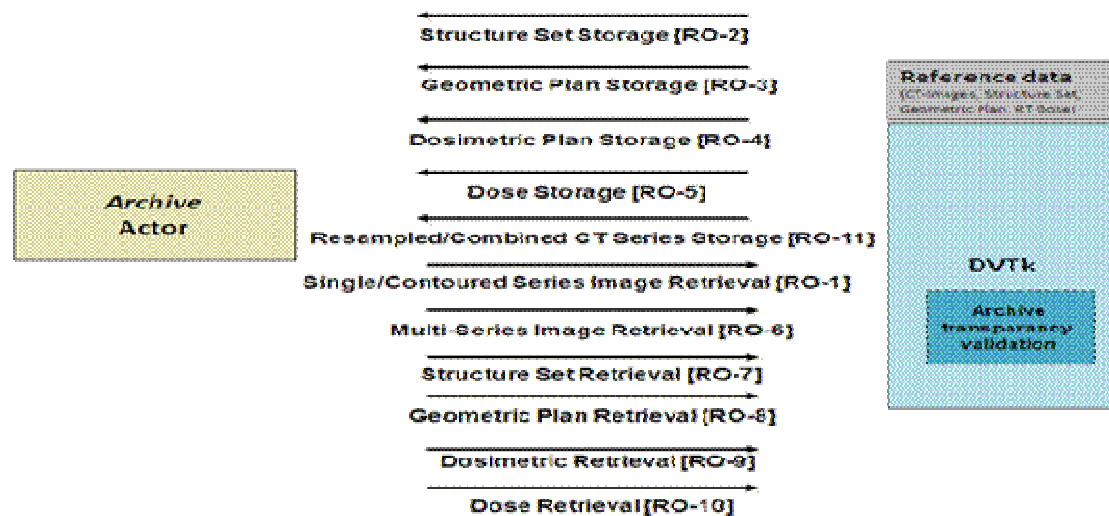

Transactions sent to actor:

- RO-2 (Structure Set Storage)
- RO-3 (Geometric Plan Storage)
- RO-4 (Dosimetric Plan Storage)
- RO-5 (Dose Storage)
- RO-11 (Resampled/Combined CT Series Storage) *

Expected transactions from actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)
- RO-8 (Geometric Plan Retrieve)
- RO-9 (Dosimetric Plan Retrieve)
- RO-10 (Dose Retrieve)

Used scenario dataset:

This scenario uses the reference test dataset.


* The single CT-Image series is used as Resampled/Combined CT Series

Test Scenario 1 Contourer actor

Description:

This test scenario provides the contourer actor with RO-1 and receives RO-2

- Checks the correct storage of Single Series CT dataset (DICOM association & interaction level)
- Validation of the created Structure set (RO-2) (including RO critical modules)
- Verifies the critical attribute mapping requirements

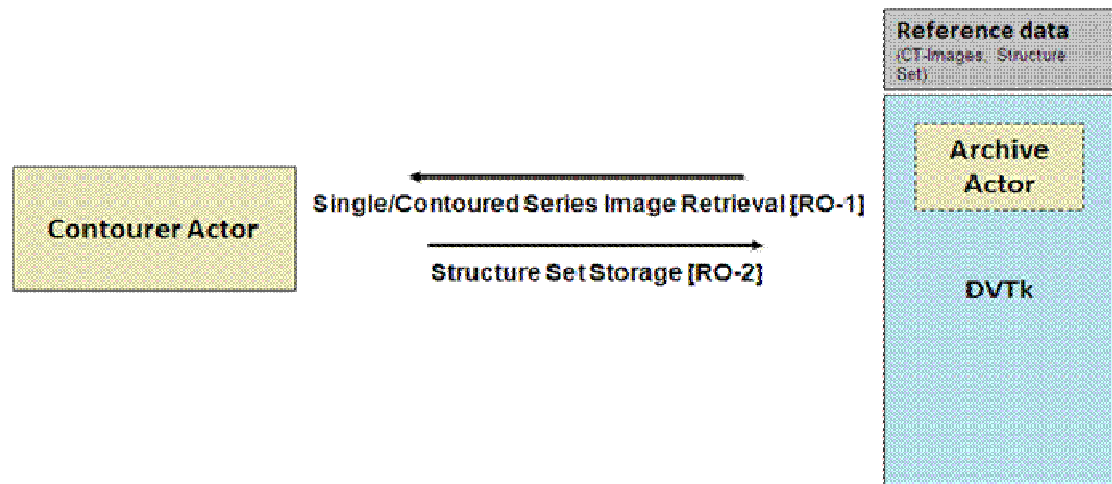

Transactions sent to actor:

RO-1 (Single/Contoured Series Image Retrieval)

Expected transactions from actor:

RO-2 (Structure Set Storage)


Used scenario dataset:

This scenario uses the reference test dataset.

Test Scenario 2 Contourer actor

Description:

- This test scenario provides the contourer actor with RO-6 and receives RO-11 and RO-2
- Checks the correct storage of Multi Series CT dataset (DICOM association & interaction level)
 - Validation of the created Structure set (RO-2) (including RO critical modules)
 - Validation of the resampled CT dataset (RO-11)
 - Verifies the critical attribute mapping requirements

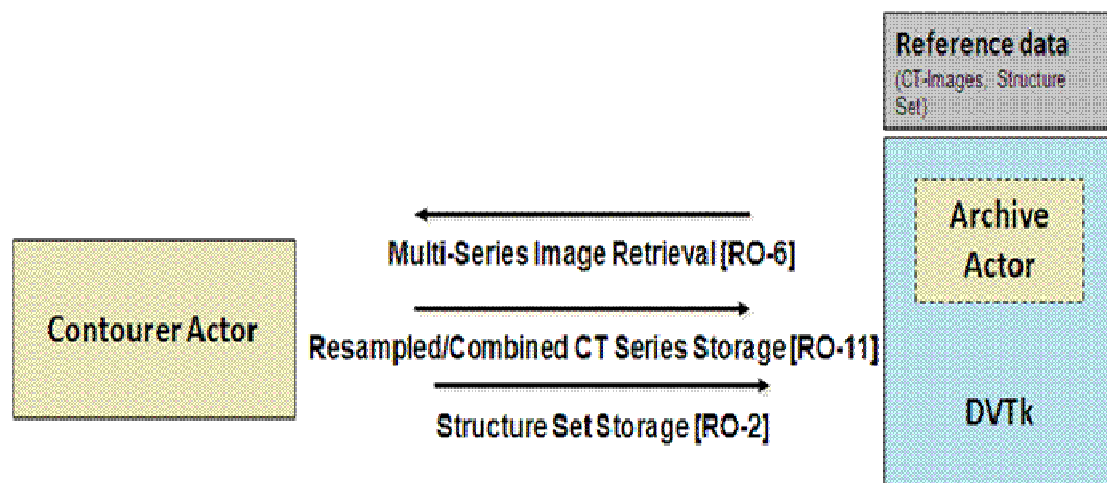

Transactions sent to actor:

RO-6 (Multi Series Image Retrieval)

Expected transactions from actor:

RO-11 (Resampled/Combined CT Series Storage)
RO-2 (Structure Set Storage)

Used scenario dataset:


This scenario uses the reference test dataset.

Test Scenario 3 Contourer actor

Description:

This test scenario provides the contourer actor with RO-1 and RO-7 and receives RO-2

- Checks the correct storage of Single Series CT dataset (DICOM association & interaction level)
- Validation of the created Structure set (RO-2) (including RO critical modules)
- Verifies the critical attribute mapping requirements

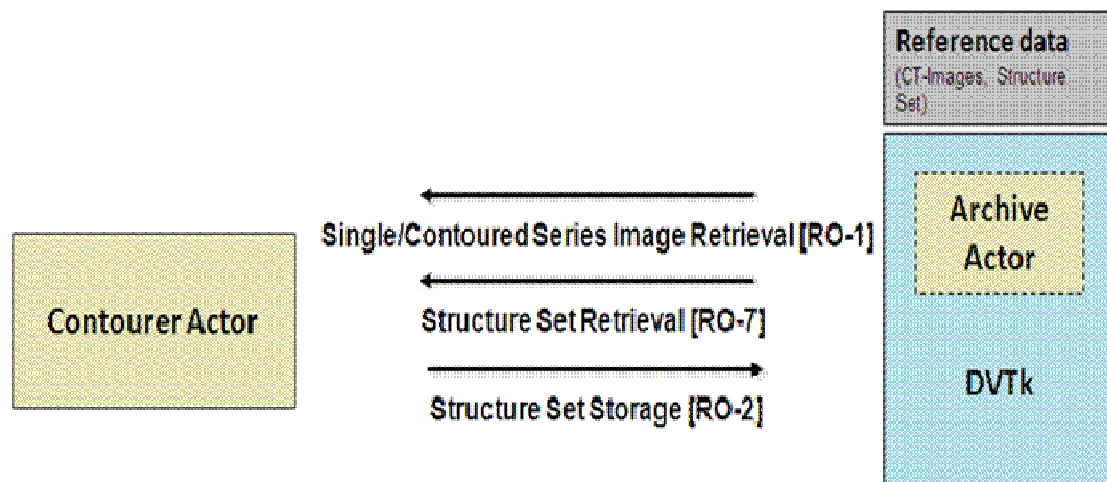

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)

Expected transactions from actor:

- RO-2 (Structure Set Storage)

Used scenario dataset:


This scenario uses the reference test dataset.

Equidistant Scenario 1 Contourer actor

Description:

This test scenario provides the contourer actor with RO-1 and receives RO-2

- Checks the correct storage of Single Series CT dataset (DICOM association & interaction level)
- Validation of the created Structure set (RO-2) (including RO critical modules)
- Verifies the critical attribute mapping requirements

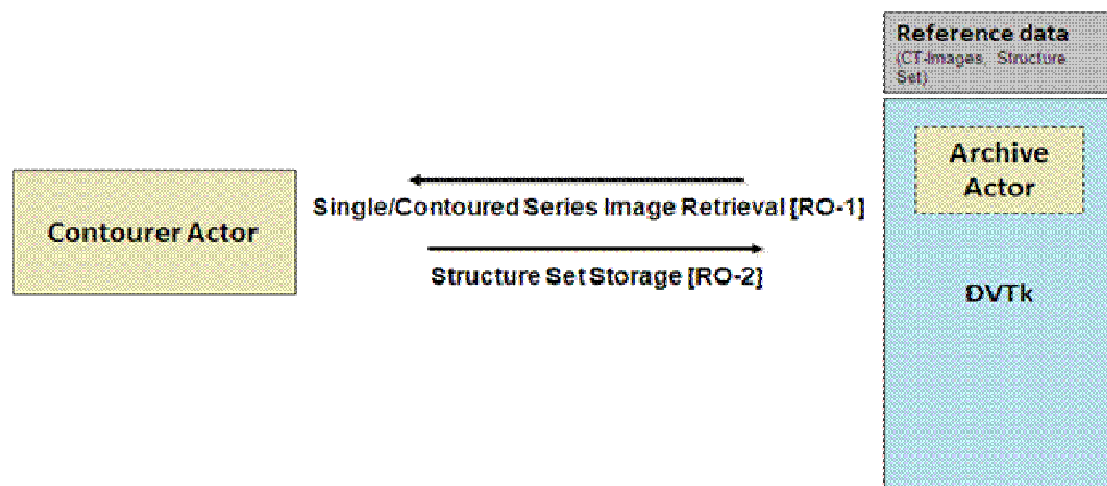

Transactions sent to actor:

RO-1 (Single/Contoured Series Image Retrieval)

Expected transactions from actor:

RO-2 (Structure Set Storage)

Used scenario dataset:

This scenario uses the equidistant reference test dataset.


Note: This test data is not clinically correct.

Test Scenario 1 Dose Displayer

Description:

This test scenario provides the dose displayer with RO-1, RO-7, RO-9 and RO-10

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Checks the correct storage of the (Dosimetric) Plan (DICOM association & interaction level)
- Checks the correct storage of RT Dose (DICOM association & interaction level)

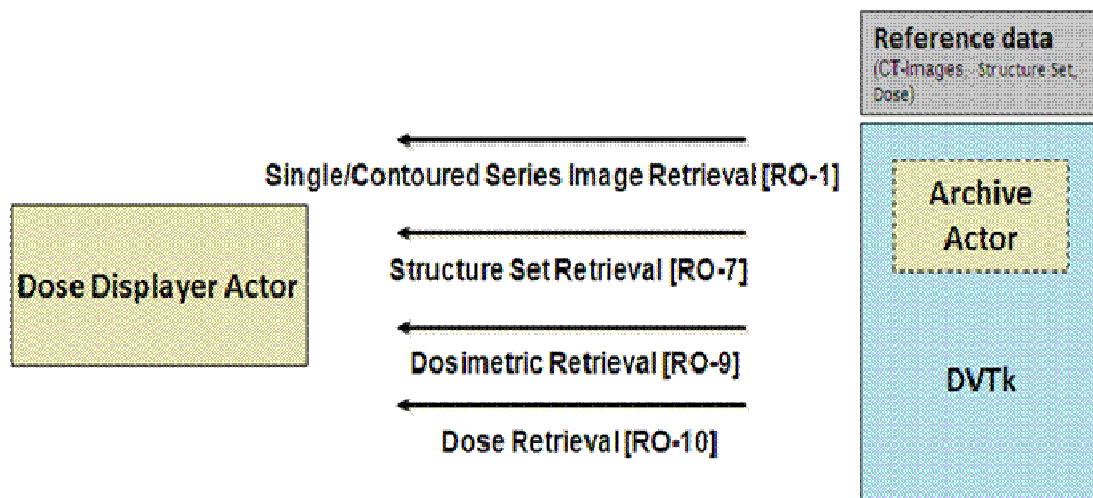

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)
- RO-9 (Dosimetric Retrieval)
- RO-10 (Dose Retrieval)

Expected transactions from actor:

n.a.

Used scenario dataset:


This scenario uses the reference test dataset.

Test Scenario 2 Dose Displayer

Description:

This test scenario provides the dose displayer with RO-1, RO-7, RO-9 and RO-10 (with minimal attributes in reference dataset (NO U modules and optional attributes))

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Checks the correct storage of the (Dosimetric) Plan (DICOM association & interaction level)
- Checks the correct storage of RT Dose (DICOM association & interaction level)

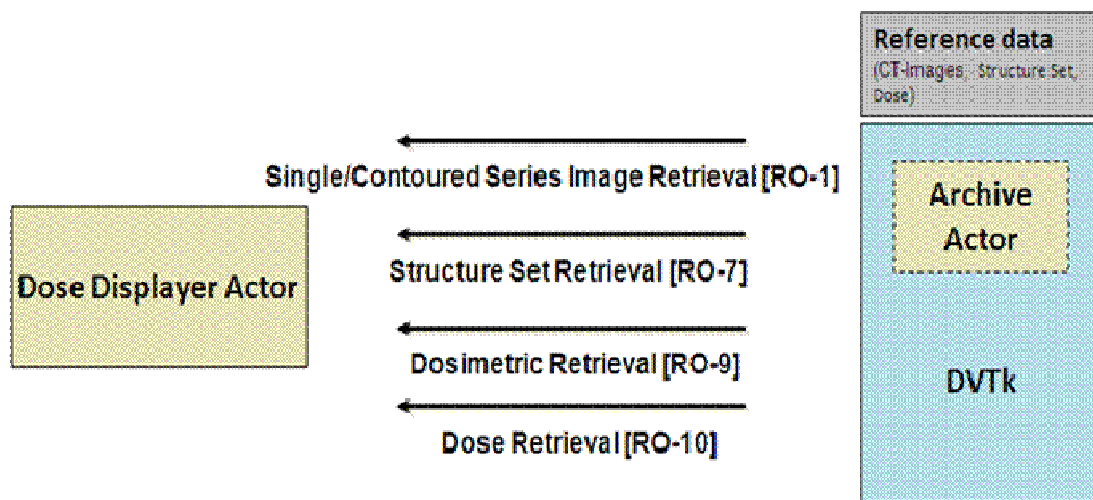

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)
- RO-9 (Dosimetric Retrieval)
- RO-10 (Dose Retrieval)

Expected transactions from actor:

n.a.

Used scenario dataset:


This scenario uses the minimal reference test dataset.

Equidistant Test Scenario 1 Dose Displayer

Description:

This test scenario provides the dose displayer with RO-1, RO-7, RO-9 and RO-10

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Checks the correct storage of the (Dosimetric) Plan (DICOM association & interaction level)
- Checks the correct storage of RT Dose (DICOM association & interaction level)

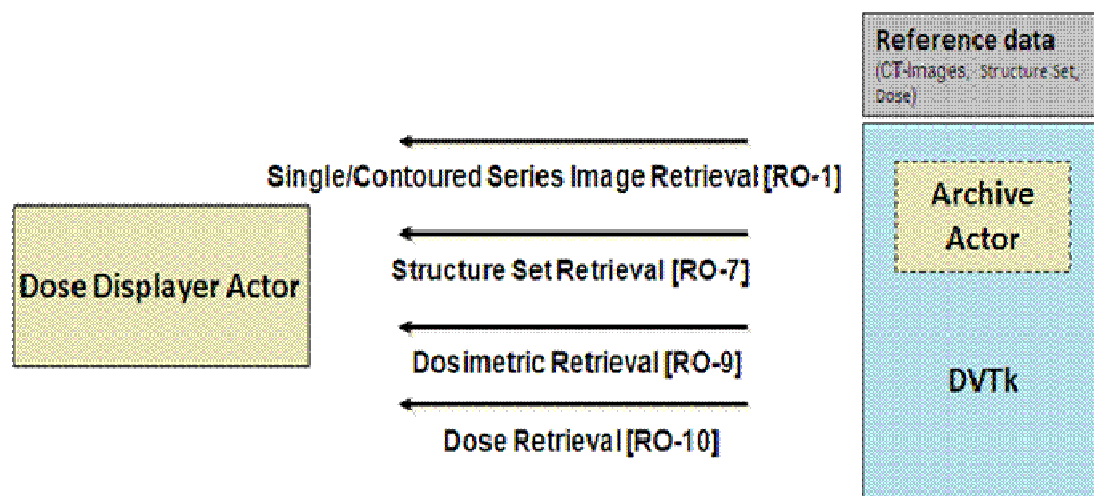

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)
- RO-9 (Dosimetric Retrieval)
- RO-10 (Dose Retrieval)

Expected transactions from actor:

n.a.

Used scenario dataset:

This scenario uses the equidistant reference test dataset.


Note: This test data is not clinically correct.

Test Scenario 1 Dosimetric Planner

Description:

This test scenario provides the dosimetric planner with RO-1, RO-7 and RO-8 and receives RO-4 and RO-5

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Checks the correct storage of the (Geometric) Plan (DICOM association & interaction level)
- Validation of the created Dosimetric Plan (RO-4) (including RO critical modules)
- Verifies the critical attribute mapping requirements

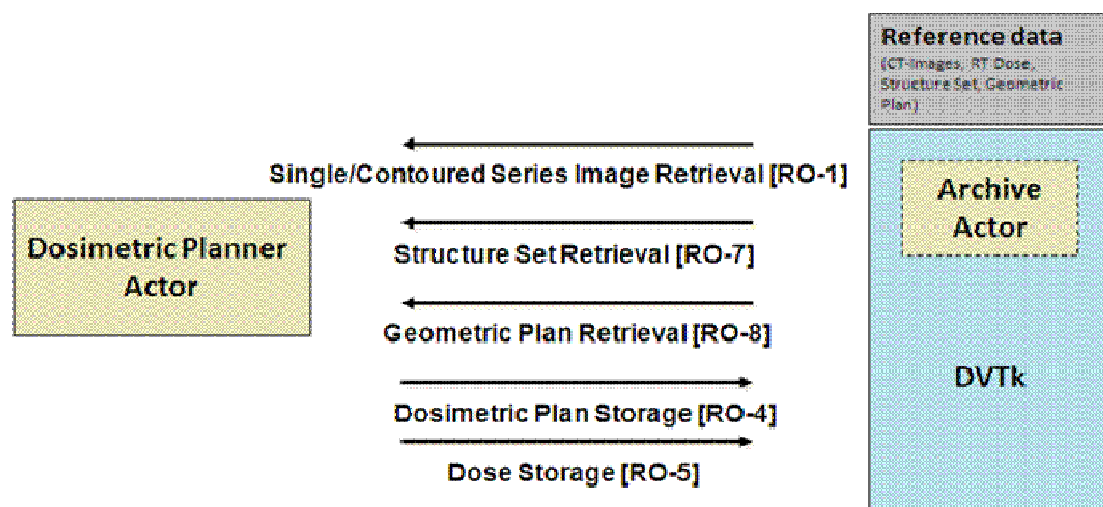

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)
- RO-8 (Geometric Plan Retrieval)

Expected transactions from actor:

- RO-4 (Dosimetric Plan Storage)
- RO-5 (Dose Storage)

Used scenario dataset:


This scenario uses the reference test dataset.

Test Scenario 2 Dosimetric Planner

Description:

This test scenario provides the dosimetric planner with RO-1, RO-7 and RO-8 and receives RO-4 and RO-5 (with minimal attributes in reference dataset (NO U modules and optional attributes))

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Checks the correct storage of the (Geometric) Plan (DICOM association & interaction level)
- Validation of the created Dosimetric Plan (RO-4) (including RO critical modules)
- Verifies the critical attribute mapping requirements

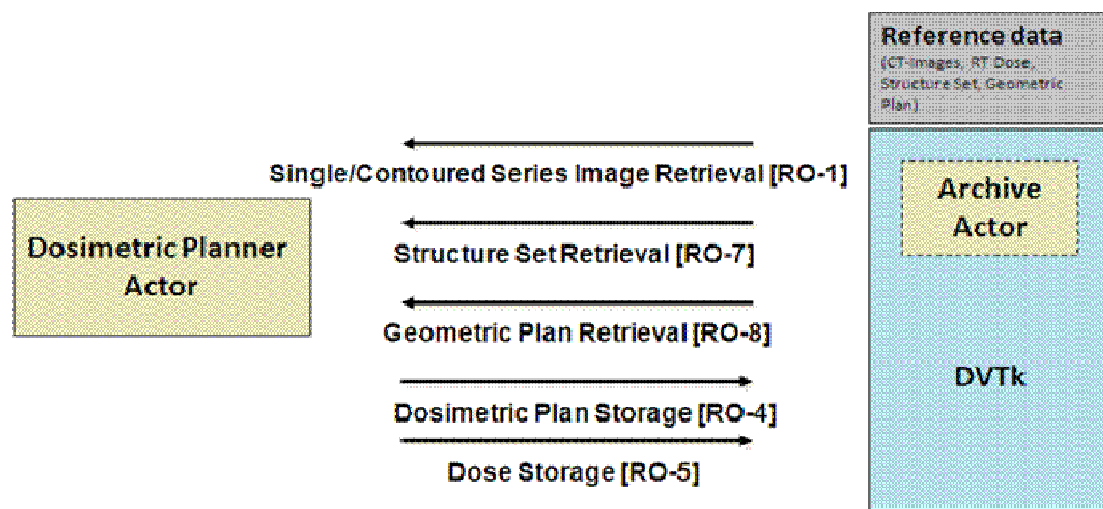

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)
- RO-8 (Geometric Plan Retrieval)

Expected transactions from actor:

- RO-4 (Dosimetric Plan Storage)
- RO-5 (Dose Storage)

Used scenario dataset:


This scenario uses the minimal reference test dataset.

Equidistant Test Scenario 1 Dosimetric Planner

Description:

This test scenario provides the dosimetric planner with RO-1, RO-7 and RO-8 and receives RO-4 and RO-5

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Checks the correct storage of the (Geometric) Plan (DICOM association & interaction level)
- Validation of the created Dosimetric Plan (RO-4) (including RO critical modules)
- Verifies the critical attribute mapping requirements

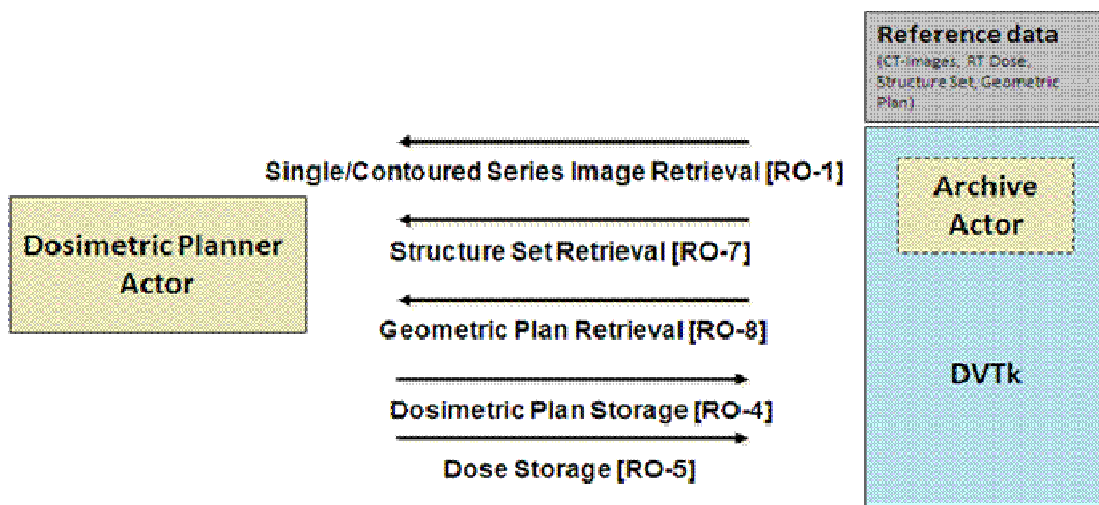

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)
- RO-8 (Geometric Plan Retrieval)

Expected transactions from actor:

- RO-4 (Dosimetric Plan Storage)
- RO-5 (Dose Storage)

Used scenario dataset:


This scenario uses the equidistant reference test dataset.
Note: This test data is not clinically correct.

Test Scenario 1 Geometric Planner

Description:

This test scenario provides the geometric planner with RO-1 and RO-7 and receives RO-3

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Validation of the created Geometric Plan (RO-3) (including RO critical modules)
- Verifies the critical attribute mapping requirements

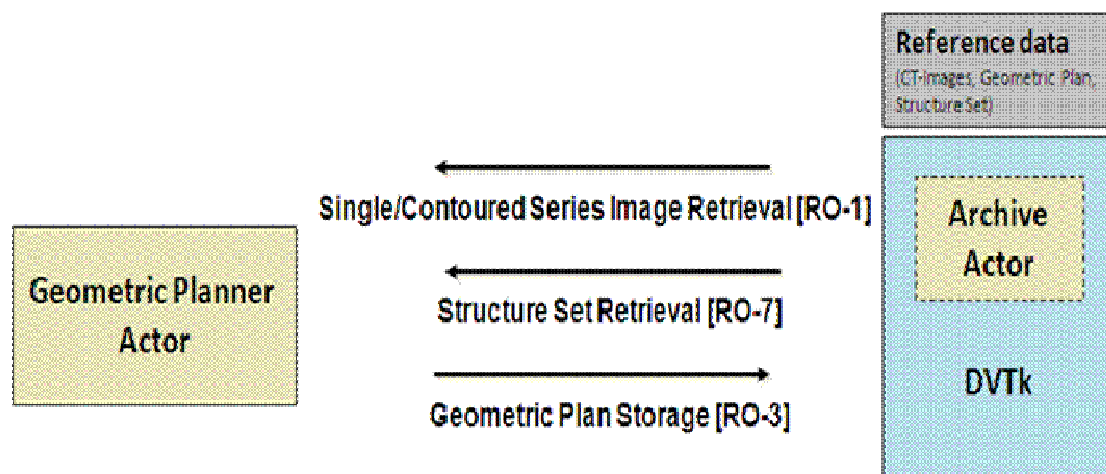

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)

Expected transactions from actor:

- RO-3 (Geometric Plan Storage)

Used scenario dataset:


This scenario uses the reference test dataset.

Test Scenario 2 Geometric Planner

Description:

This test scenario provides the geometric planner with RO-1 and RO-7 and receives RO-3 (with minimal attributes in reference dataset (NO U modules and optional attributes))

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Validation of the created Geometric Plan (RO-3) (including RO critical modules)
- Verifies the critical attribute mapping requirements

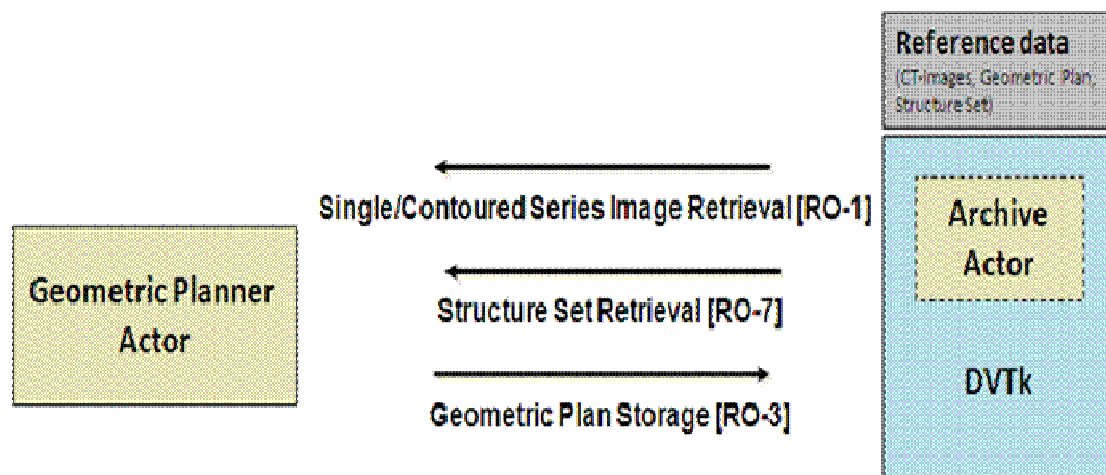

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)

Expected transactions from actor:

- RO-3 (Geometric Plan Storage)

Used scenario dataset:


This scenario uses the minimal reference test dataset.

Equidistant Test Scenario 1 Geometric Planner

Description:

This test scenario provides the geometric planner with RO-1 and RO-7 and receives RO-3

- Checks the correct storage of Single/contoured CT series dataset (DICOM assoc. & interaction level)
- Checks the correct storage of the Structure set (DICOM association & interaction level)
- Validation of the created Geometric Plan (RO-3) (including RO critical modules)
- Verifies the critical attribute mapping requirements

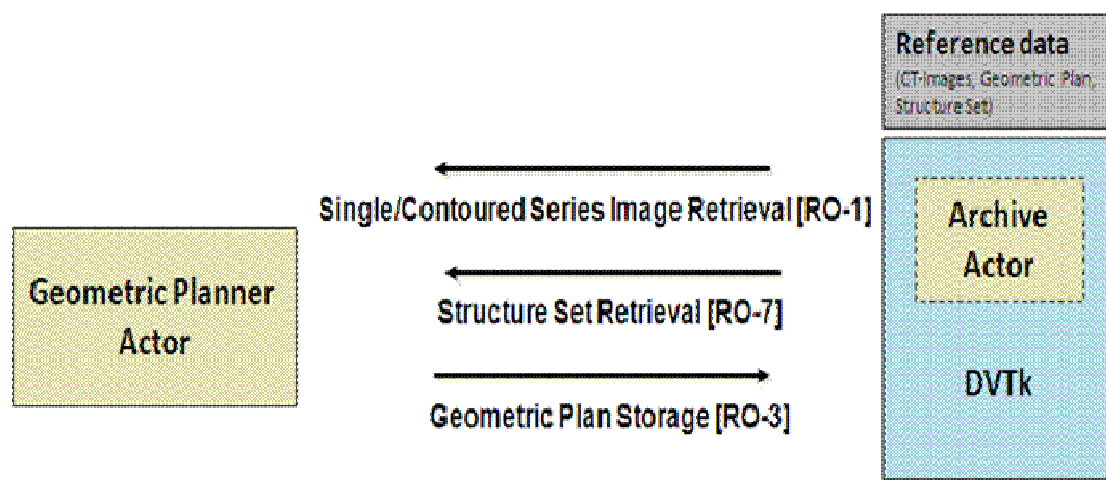

Transactions sent to actor:

- RO-1 (Single/Contoured Series Image Retrieval)
- RO-7 (Structure Set Retrieval)

Expected transactions from actor:

- RO-3 (Geometric Plan Storage)

Used scenario dataset:

This scenario uses the equidistant reference test dataset.

Note: This test data is not clinically correct.